

Erasmus NC

0.018" system

Non-Compliant OTW

Excellent for high pressure dilatation and lesion preparation

Combination of strength and optimal control in routine and challenging cases

Broad range of sizes available

Confident to treat all lower extremity lesions

Tri-fold (3 pleats) for \varnothing 2.0-4.0mm

Five-fold (5 pleats) for \varnothing 4.5-8.0mm

Seamless tip-to-shaft design


Provides smooth transmission and maximized pushability

Low tip profile 0.022"

Facilitates lesion entry

Optimized non-compliant balloon material

High pressure balloon inflation up to 22atm

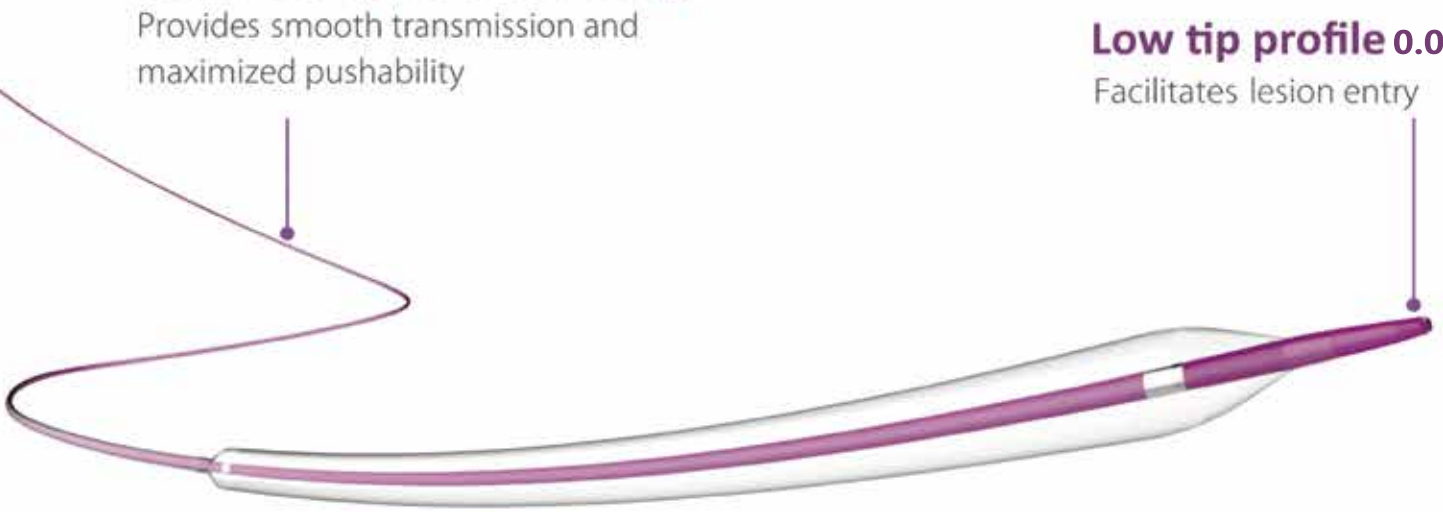

Specifications

| Balloon Diameter (mm) | Balloon Length (mm) | | | | | | | | | Catheter Working Length (cm) | Sheath Compatibility (F) | Nominal Pressure (atm) | Rated Burst Pressure (atm) |
|-----------------------|---------------------|----|----|----|----|----|----|-----|-----|------------------------------|--------------------------|------------------------|----------------------------|
| | 10 | 15 | 20 | 30 | 40 | 60 | 80 | 100 | 150 | | | | |
| 3.0 | | | V | V | V | V | | | | 40 | 4 | 12 | 22 |
| 3.5 | | | V | V | V | V | | | | 40 | 4 | 12 | 22 |
| 4.0 | | | V | V | V | V | | | | 40 | 4 | 12 | 22 |
| 4.5 | | | V | V | V | V | | | | 40 | 5 | 12 | 20 |
| 5.0 | | | V | V | V | V | | | | 40 | 5 | 12 | 20 |
| 5.5 | | | V | V | V | V | | | | 40 | 5 | 12 | 20 |
| 6.0 | | | V | V | V | V | | | | 40 | 6 | 12 | 20 |
| 7.0 | | | V | V | V | V | | | | 40 | 6 | 12 | 18 |
| 8.0 | | | V | V | V | V | | | | 40 | 6 | 12 | 18 |
| 2.0 | V | V | V | V | V | V | V | V | V | 70 | 4 | 12 | 22 |
| 2.5 | V | V | V | V | V | V | V | V | V | 70 | 4 | 12 | 22 |
| 3.0 | V | V | V | V | V | V | V | V | V | 70 | 4 | 12 | 22 |
| 3.5 | V | V | V | V | V | V | V | V | V | 70 | 4 | 12 | 22 |
| 4.0 | V | V | V | V | V | V | V | V | V | 70 | 4 | 12 | 22 |
| 4.5 | V | V | V | V | V | V | V | V | V | 70 | 5 | 12 | 20 |
| 5.0 | V | V | V | V | V | V | V | V | V | 70 | 5 | 12 | 20 |
| 5.5 | V | V | V | V | V | V | V | V | V | 70 | 5 | 12 | 20 |
| 6.0 | V | V | V | V | V | V | V | V | V | 70 | 6 | 12 | 20 |
| 7.0 | V | V | V | V | V | V | V | V | V | 70 | 6 | 12 | 18 |
| 8.0 | V | V | V | V | V | V | V | V | V | 70 | 6 | 12 | 18 |
| 2.0 | V | V | V | V | V | V | V | V | V | 90 | 4 | 12 | 22 |
| 2.5 | V | V | V | V | V | V | V | V | V | 90 | 4 | 12 | 22 |
| 3.0 | V | V | V | V | V | V | V | V | V | 90 | 4 | 12 | 22 |
| 3.5 | V | V | V | V | V | V | V | V | V | 90 | 4 | 12 | 22 |
| 4.0 | V | V | V | V | V | V | V | V | V | 90 | 4 | 12 | 22 |
| 4.5 | V | V | V | V | V | V | V | V | V | 90 | 5 | 12 | 20 |
| 5.0 | V | V | V | V | V | V | V | V | V | 90 | 5 | 12 | 20 |
| 5.5 | V | V | V | V | V | V | V | V | V | 90 | 5 | 12 | 20 |
| 6.0 | V | V | V | V | V | V | V | V | V | 90 | 6 | 12 | 20 |
| 7.0 | V | V | V | V | V | V | V | V | V | 90 | 6 | 12 | 18 |
| 8.0 | V | V | V | V | V | V | V | V | V | 90 | 6 | 12 | 18 |
| 2.0 | V | V | V | V | V | V | V | V | V | 150 | 4 | 12 | 22 |
| 2.5 | V | V | V | V | V | V | V | V | V | 150 | 4 | 12 | 22 |
| 3.0 | V | V | V | V | V | V | V | V | V | 150 | 4 | 12 | 22 |
| 3.5 | V | V | V | V | V | V | V | V | V | 150 | 4 | 12 | 22 |
| 4.0 | V | V | V | V | V | V | V | V | V | 150 | 4 | 12 | 22 |
| 4.5 | V | V | V | V | V | V | V | V | V | 150 | 5 | 12 | 20 |
| 5.0 | V | V | V | V | V | V | V | V | V | 150 | 5 | 12 | 20 |
| 5.5 | V | V | V | V | V | V | V | V | V | 150 | 5 | 12 | 20 |
| 6.0 | V | V | V | V | V | V | V | V | V | 150 | 6 | 12 | 20 |
| 7.0 | V | V | V | V | V | V | V | V | V | 150 | 6 | 12 | 18 |
| 8.0 | V | V | V | V | V | V | V | V | V | 150 | 6 | 12 | 18 |

Remark: Sizes within have longer delivery time.



cNovate Medical BV
 Address: Hogebrinkerweg 31, 3871 KM Hoevelaken, The Netherlands
 Tel: +31 850140404 E-mail: cs@cNovate.eu Web: www.cNovate.eu

